


INTERNAL **DAS MU**




INTERNAL MULTI-UNIT DAS SYSTEM



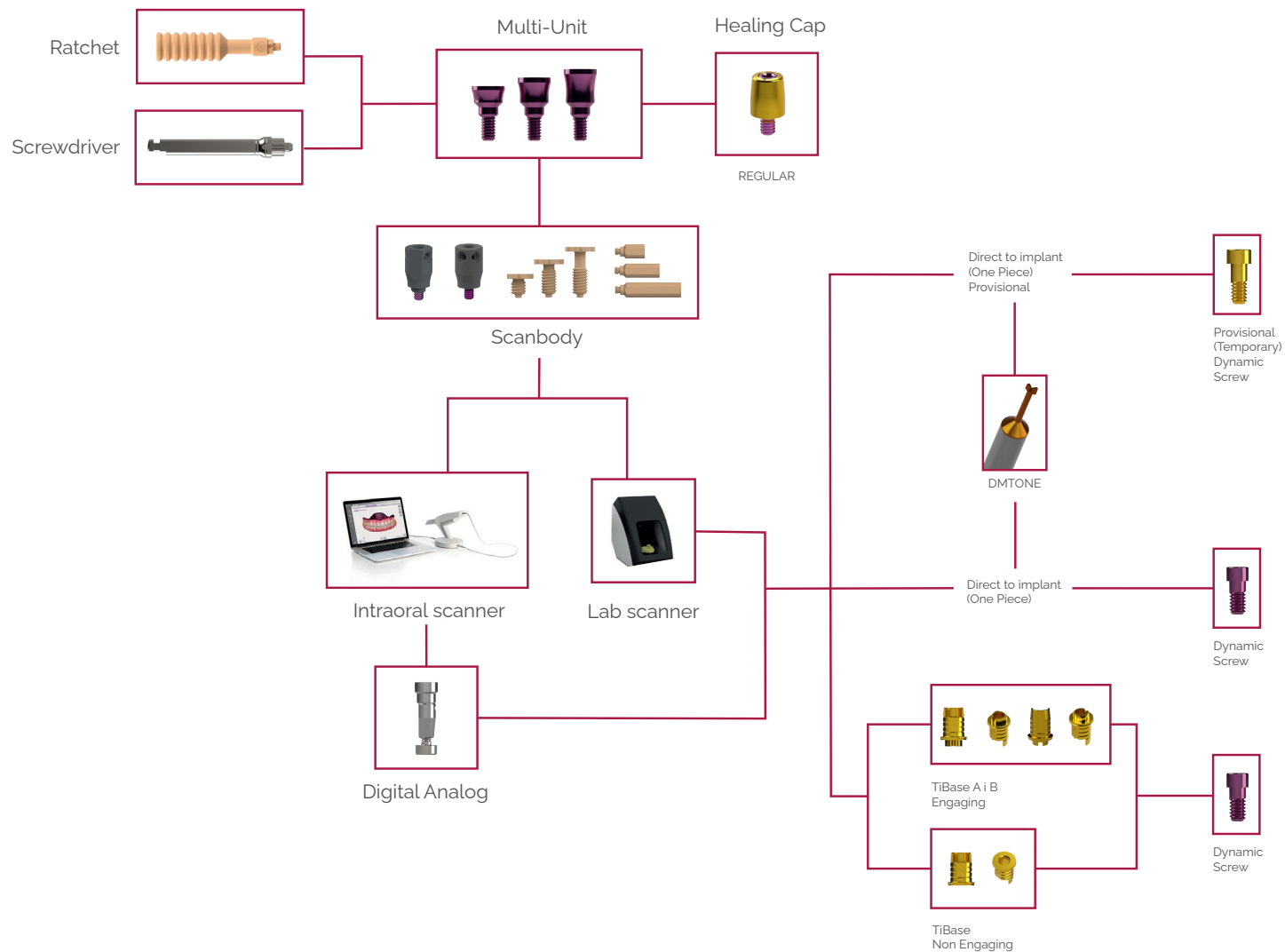

MULTI-UNIT WORKFLOW

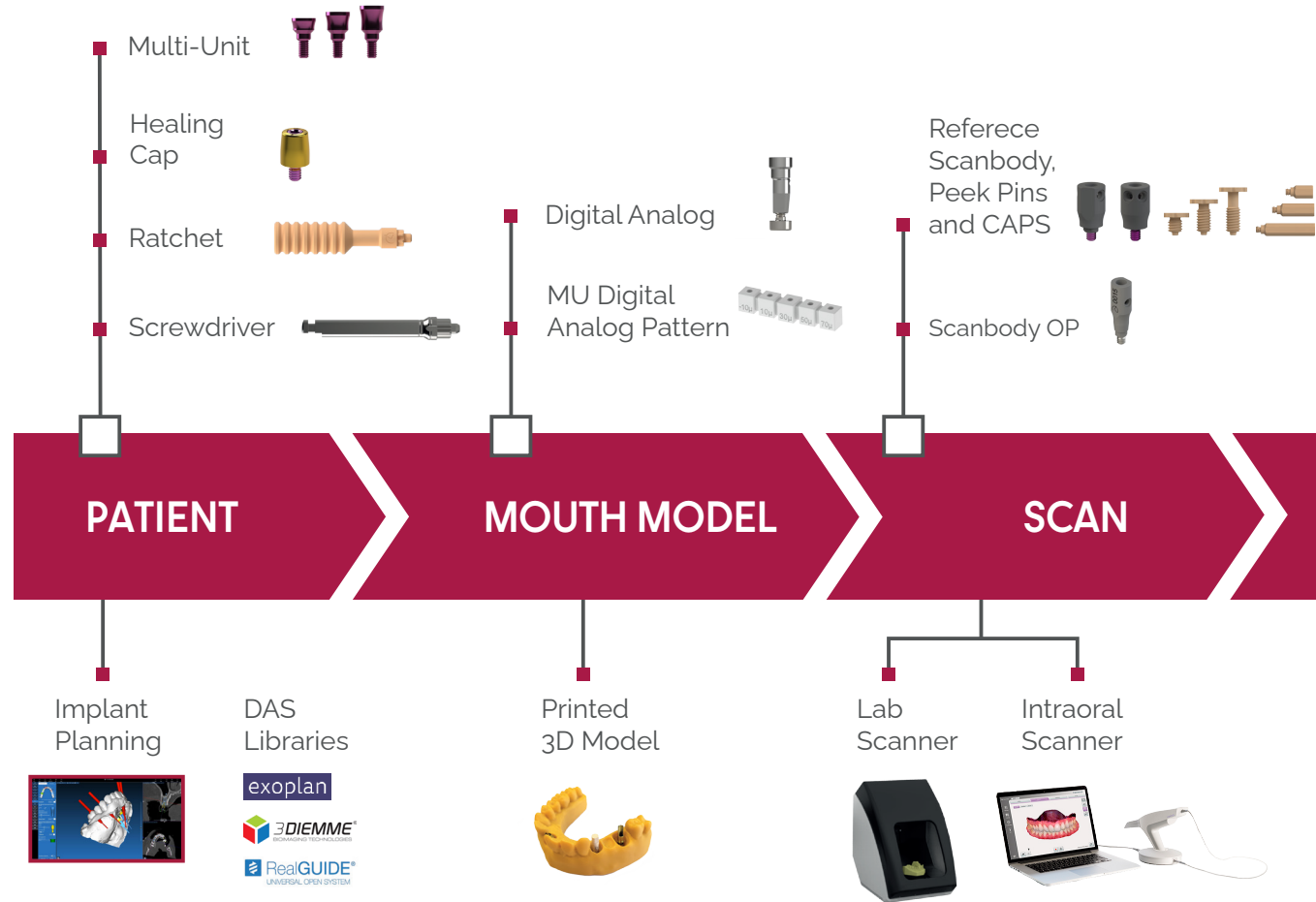

Dynamic
TiBase A i B
Engaging



TiBase
Non Engaging



Direct to Implant
(One Piece)
Provisional



Direct to Implant
(One Piece)



DMTONE



CAD Services



CAM Strategies



MU TiBase
Pattern



Dynamic Screwdriver
Dynamic Screw
Provisional Dynamic Screw



Dynamic TiBase



Direct to Implant
(One Piece)
Provisional



Direct to Implant
(One Piece)
Final Structure



DESIGN

MILLING

FINAL STRUCTURE

CAD
Design



DAS
Libraries

3shape
exocad



Milling
CAM



DAS Tested
CAM Software

worknc
Dental
hyperDENT

FOLLOW-ME!
TECHNOLOGY GROUP




GO2dental

vhf


REFERENCE SCANBODY & PEEK PINS

The Reference Scanbody is made of titanium and is designed to improve scanning for **edentulous cases**.

- Peek Pins can be used for **best scanning**
- The Reference Scanbody is compatible with different **Multi-Unit**

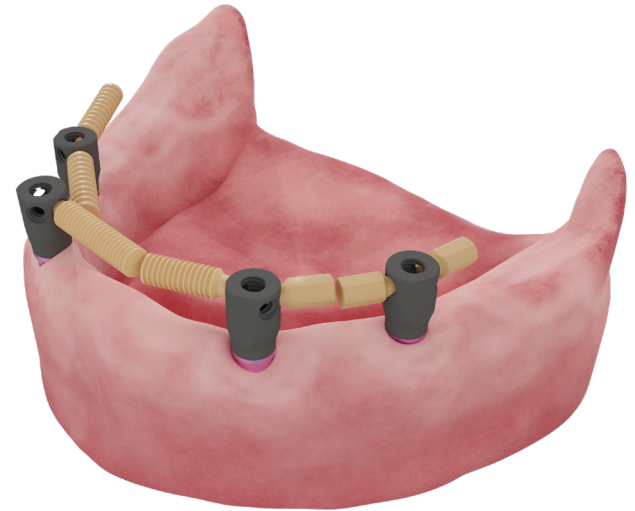

Peek Pins lengths

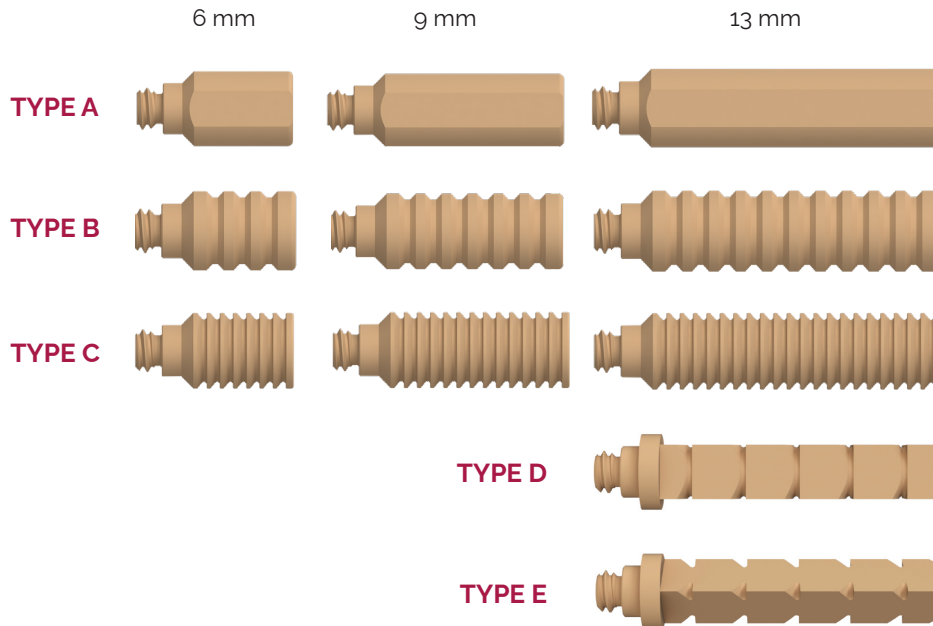

Reference Scanbody



CAPS COMPLETE ARCH PILLAR SYSTEM

The Complete Arch Pillar System (CAPS) is made to improve **scan accuracy** and to **record maxilomandibular relationship** for the fabrication of a complete arch.



With the CAPS System we can finally close the Circle of **Digital Workflow** for a Complete Arch.

CAPS heights

REGULAR CAPS



WIDE CAPS

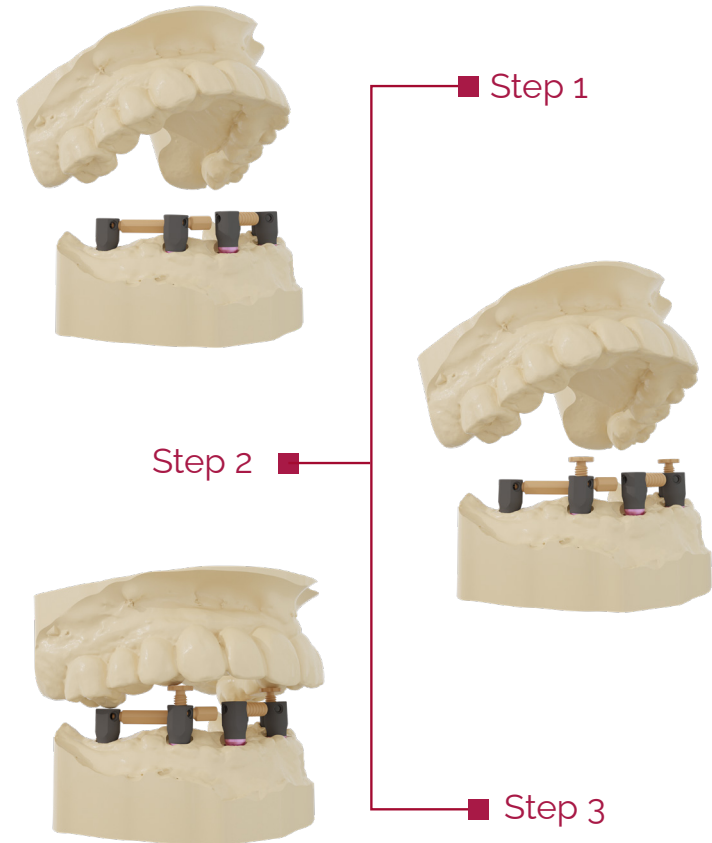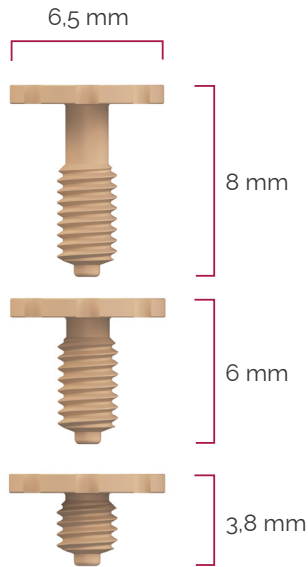

INTERNAL MULTI-UNIT

The Multi-Unit (MU) abutment has been carefully designed to rehabilitate partially or fully edentulous arches, as well as individual.

Available in different gingival heights



Maximum diameter of the MU is 4 mm.
By being narrower it avoids contact with the bone



Similar color to the gingiva
for perfect adaptation/
integration



Concave design facilitates healing
and soft tissue adaptation

**Divergence between
implants: 110°**

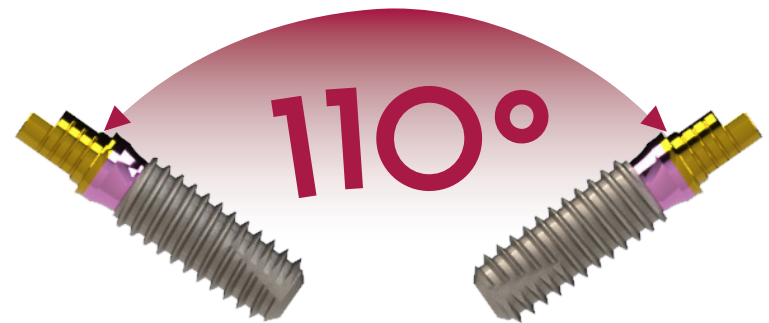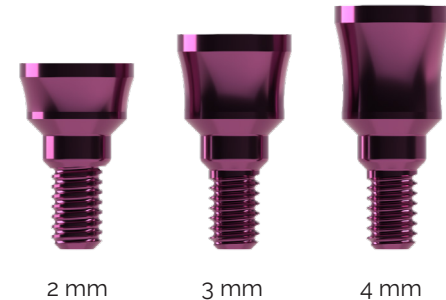

DAS MU & INTERNAL DAS MU

Using DAS MU in combination
with INTERNAL DAS MU

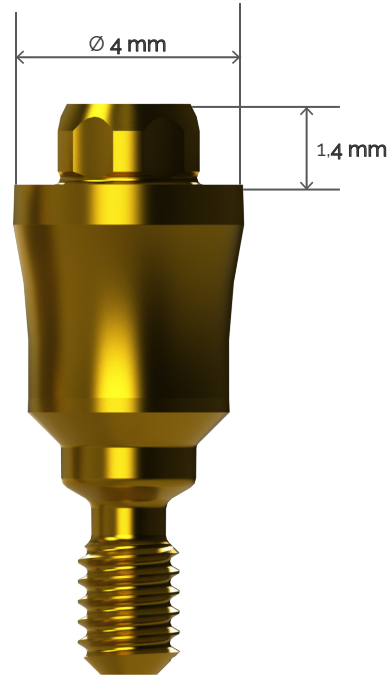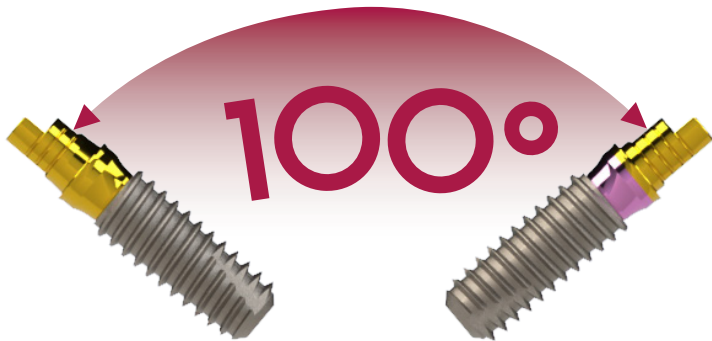

DAS MU


1 mm 2 mm 3 mm 4 mm




INTERNAL DAS MU

2 mm 3 mm 4 mm

MU INSERTION PROCEDURE

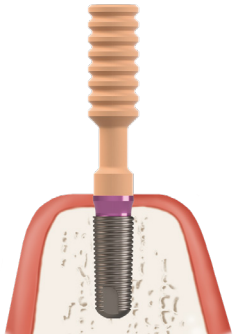

STEP 1:
Use the plastic gripping tool to attach the MU to the implant

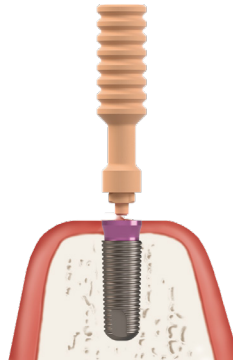


STEP 2:
Remove the plastic tool




STEP 3:
Tighten MU to **30 N.cm** using the MU Screwdriver

MU HEALING CAPS

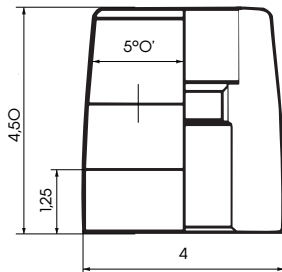

HEALING CAP REGULAR

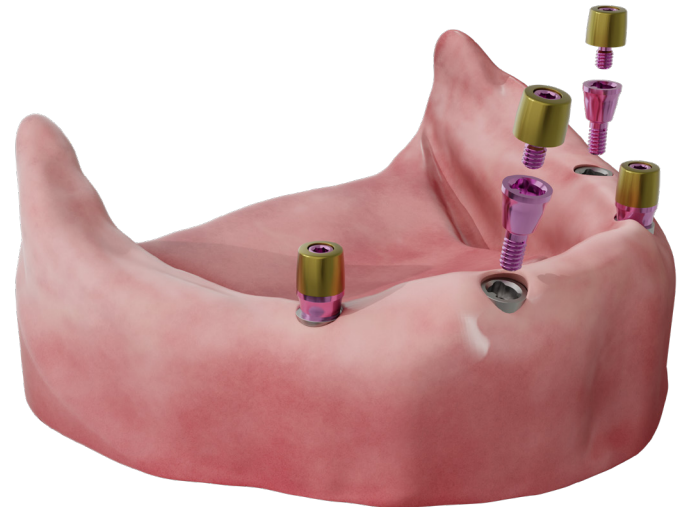


DIGITAL ANALOG

Digital analog of the dental implant to simulate implant position in a 3D printed dental model.




Concave notch

Top precision in longitudinal position

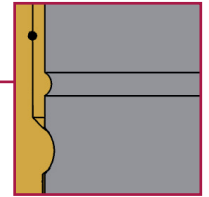


Screwed fastening

Prevents the analog from moving in Z

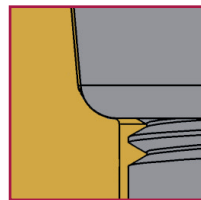

Curved surface

Accuracy of orientation guaranteed

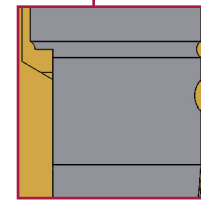

Longitudinal cut

Longitudinal cut to avoid rotation X-Y

DYNAMIC TIBASE

TO CORRECT ANGULATION



Connection slot

The window is oriented
at 0° with respect to
the slot

Connection slot

The window is oriented
at 180° with respect to
the slot


TIBASE
TYPE A
NR=ENGAGING

TIBASE
TYPE B
NR=ENGAGING

TIBASE
R=NON ENGAGING

Angulated Screw
Channel

CAD DESIGN

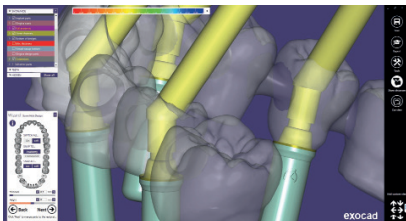


3shape
exocad
dental wings

Internal DAS Multi-Unit system allows to change the MU in the library **without rescanning or redesigning** the case. That eases the lab and clinic work, as technician can change the MU without the need to make a new appointment with the patient to re-scan.


If we make a NR case, the option of changing MU is not allowed in the design.

DIRECT TO IMPLANT ONE PIECE

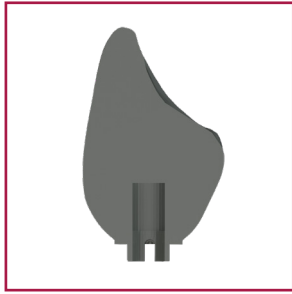


Step 1

Crown with pre drill

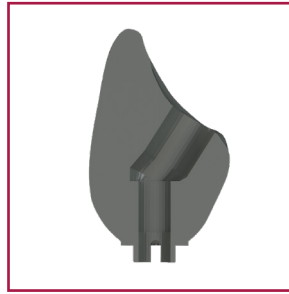



Step 2

Crown with angled channel

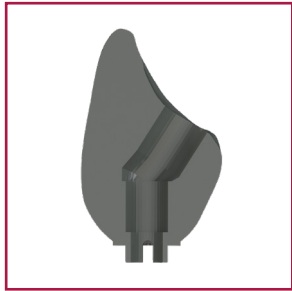

Step 3 DMTone

Crown with Dynamic Milling Tool. Milling the screw seat and increasing the diameter of the straight channel

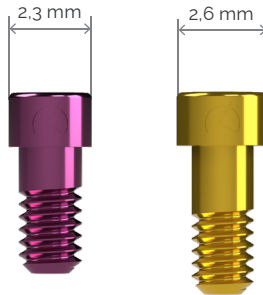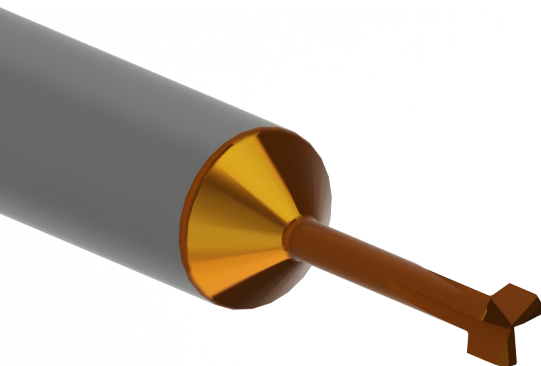

Dynamic screw

Provisional (temporary) Dynamic screw

DMTONE
DYNAMIC MILLING TOOL



Final restorations 0-35°

Direct to implant libraries available R/NR

R = Non engaging | NR = Engaging



Greater seating thickness

Greater contact surface between the screw and the material for greater safety in soft materials

Provisional (temporary) restorations up to 35°

Provisional (temporary) direct to implant libraries available R/NR

R = Non engaging | NR = Engaging

INTERNAL DAS MU SYSTEM COMPONENTS



Ratchet
49.409.000.02-2



Screwdriver
43.321.320.01-2



Digital Analog
34.312.210.01-2




Digital Analog Positioner
49.309.000.01-2




Dynamic Screwdriver
43.618.201.01-2 (18 mm)
43.624.201.01-2 (24 mm)
43.632.201.01-2 (32 mm)




MU DMTone
33.390.716.01-2 Shank Ø3
33.490.716.01-2 Shank Ø4
33.690.716.01-2 Shank Ø6

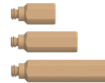

Peek Pins
TYPE A
49.414.000.01-2 (6 mm)
49.415.000.01-2 (9 mm)
49.416.000.01-2 (13 mm)

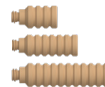

TYPE B
49.414.000.02-2 (6 mm)
49.415.000.02-2 (9 mm)
49.416.000.02-2 (13 mm)

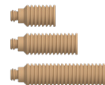

TYPE C
49.414.000.03-2 (6 mm)
49.415.000.03-2 (9 mm)
49.416.000.03-2 (13 mm)



TYPE D
49.416.000.04-2 (13 mm)



TYPE E
49.416.000.05-2 (13 mm)



CAPS
Regular
49.418.000.01-2 (3,8 mm)
49.419.000.01-2 (6 mm)
49.420.000.01-2 (8 mm)



Wide
49.418.000.02-2 (3,8 mm)
49.419.000.02-2 (6 mm)
49.420.000.02-2 (8 mm)



Reference Scanbody
Engaging
54.312.210.31-2



Non-engaging
54.322.210.31-2



Healing Cap
Regular
40.318.003.89-2.



MU Dynamic TiBase
Type A Engaging
31.312.210.01-2



Type B Engaging
31.313.210.01-2



Non-engaging
31.322.210.01-2



Dynamic Screw
41.318.050.01-2



Provisional (temporary) Dynamic Screw
41.318.057.01-2

TALLADIUM GUARANTEE

TERMS AND CONDITIONS

These guarantee terms and conditions ("T&C") cover the entire range of Talladium products ("Products"), manufactured by TALLADIUM ESPAÑA S.L. and distributed by Geoda Medical S.L. or official dealers. The guarantee described in these T&C is exclusively in benefit of the clinician ("Clinician") and of the dental technician ("Technician") and not for the benefit of third parties or institutions, including patients.

GUARANTEE PERIOD

TALLADIUM ESPAÑA S.L. offers a lifelong guarantee for its entire range of products starting from the date of issue of the invoice.

GUARANTEE SCOPE

Subject to the limitations and exceptions described in these T&C, TALLADIUM ESPAÑA S.L. will offer the following benefits:

QUALITY: If there are defects in the materials or in the manufacturing of the Product, TALLADIUM ESPAÑA S.L. will replace the Product with no additional cost.

SAFETY: If, having complied with all the product indications, the prosthesis should have to be made again, due to a fault in the Dynamic Abutment or Dynamic Titanium Base system, TALLADIUM ESPAÑA S.L. will replace the abutments and screws necessary to remake the prosthesis, as well as the costs derived from its manufacturing.

In case of having used our products and having complied with all the product indications, the implants suffer any damage, TALLADIUM ESPAÑA S.L. will pay the cost of the implants. This coverage will only be valid during the first 6 months after the collocation of the prosthesis which includes our products.

CLAIM REQUIREMENTS AND PROCEDURE

To receive the benefits indicated in these T&C, the treating Clinician must satisfy the following requirements:

- The claim must be notified to TALLADIUM ESPAÑA S.L. within (30) days since the date the claimed defect was detected.
- This requires that the Clinician or Technician must contact the customer service department by telephone or by e-mail to make the claim.
- A claim form will be completed, which, together with a document or report which justifies the faulty Product and the faulty Product itself, will be sent by the customer to TALLADIUM ESPAÑA S.L. offices, within the previously indicated period.
- Clinicians or Technicians presenting a claim in agreement with these T&C must be up to date in any payments owing to TALLADIUM ESPAÑA S.L. or to any of its subsidiaries, at the time when the claim form is presented.
- All the use procedures of our Products must be carried out in agreement with the instructions of TALLADIUM ESPAÑA S.L. as well as in accordance with commonly accepted dentistry practices.
- The expenses derived from this procedure will be assumed by the customer. The return shipping costs will be assumed by TALLADIUM ESPAÑA S.L. in all those cases covered by these T&C. Regardless of the guarantee rights, claims should be notified as soon as possible in order to comply with regulatory requirements.

GENERAL LIMITATIONS OF THIS GUARANTEE

With the exception of the guarantee described in these T&C, neither TALLADIUM ESPAÑA S.L. nor its representatives, nor third parties manufacturing or distributing the Products, represent or offer a guarantee, agreement or any other express or implicit, oral or written, commitment, with respect to the Products (without limitation), including guarantees involved in the marketing, durability or suitability for individual uses or purposes. In addition and within the maximum extent permitted by the relative law, TALLADIUM ESPAÑA S.L. rejects (on its own behalf, and on behalf of its representatives and third parties that manufacture or distribute Products) any responsibility with respect to any direct or indirect damage caused, which may result from or be a consequence of the design, composition of the dental prosthesis into which the Products are integrated.


GUARANTEE EXCLUSIONS

TALLADIUM ESPAÑA S.L. limits this guarantee to:

- Transformed abutments that form part of the dental prosthesis. But not the screws used to anchor them.
- Those products that are not used with the accessories and parts marketed by Talladium España

AMENDMENT OR SUSPENSION OF THE GUARANTEE


TALLADIUM ESPAÑA S.L. reserves the right to amend or withdraw these T&C at any time and without prior notification. Any modification or suspension shall not affect products already placed in patients.

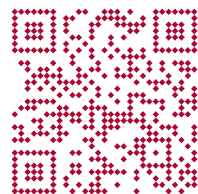

This printed version is probably
outdated and may not include all
products or compatibilities.

**Scan and download the
latest version of this catalogue**

www.dynamicabutment.com

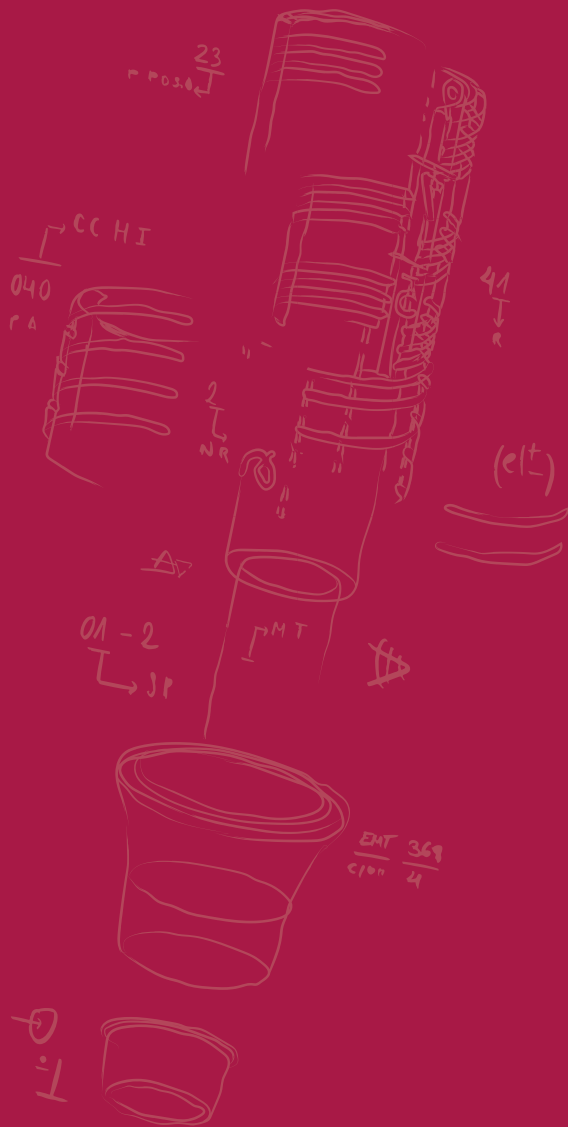

**DYNAMIC
ABUTMENT
SOLUTIONS**

Virginia Woolf, 17
25005 · Lleida (Spain)

INTERNATIONAL +34 873 450 709
das@dynamicabutment.com

SPAIN +34 973 289 580
spain@dynamicabutment.com



www.dynamicabutment.com

